

Genumedi® PT

Knee support promotes muscular balance in cases of patellar lateralisation




Indications

- patellofemoral pain syndrome in cases of:
 - hypermobile patella
 - patellar lateralisation
 - muscular imbalance
 - patellar maltracking
- post-operative treatment (patella stabilisation)


Mode of action

- positive effect on the muscular structures that guide the patella to compensate for imbalance
- reduces muscle tone of vastus lateralis
- activates vastus medialis
- laterally reinforced silicone ring gently guides the patella for reliable movement in the femoral groove
- proprioceptive compression for reliable guidance and stabilisation of the joints

Product benefits

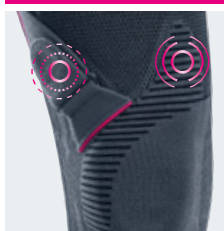

Relaxes and activates
Precisely placed 3D pads relax and activate the vastus lateralis and vastus medialis muscles respectively.

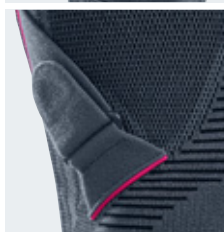

Individual adjustment
The control strap and thus the pressure on the vastus lateralis muscle can be adjusted as necessary.




Correcting patellar tilt
The silicone ring surrounding the patella is laterally reinforced to promote correct trochlear articulation. Patellar tilt is corrected.



Clima Comfort




Clima Fresh

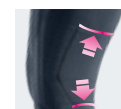

Comfort Zone



medi Kidz – in sizes I + II also for children suitable.

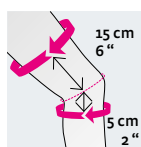

circumference "f" (cm)	37 – 40	40 – 43	43 – 46	46 – 49	49 – 52	52 – 55	55 – 58
circumference "d" (cm)	25 – 28	28 – 31	31 – 34	34 – 37	37 – 40	40 – 43	43 – 46
size	I	II	III	IV	V	VI	VII
article number	K.14◇.221/002	K.14◇.222/002	K.14◇.223/002	K.14◇.224/002	K.14◇.225/002	K.14◇.226/002	K.14◇.227/002

◇ = fill in for orientation: 2= left, 3= right
Available with or without physioband and exercise poster.