

M.4s[®] PCL dynamic

4-point rigid frame brace with extension / flexion limitation for treating injuries of the posterior cruciate ligament

Further product for posterior cruciate ligament therapy: medi PTS (p. 95)

Indications

- conservative treatment after injury of the posterior cruciate ligament
- post-operative treatment after reconstruction of the posterior cruciate ligament


Mode of action

- individually adjustable protection pad actively exerts a force directed anteriorly to hold the tibia in the desired position. This protects and relieves tension on the posterior cruciate ligament.
- reliable stabilisation of the knee joint by the 4-point rigid frame construction
 - flexion and extension limitation protect against undesirable joint movements
 - extension limitation at: 0°, 10°, 20°, 30°, 45°
 - flexion limitation at: 0°, 10°, 20°, 30°, 45°, 60°, 75°, 90°
 - immobilisation at: 0°, 10°, 20°, 30°, 45°


Product benefits



PCL protection pad
the pad presses actively on the tibia and holds it in the desired position.

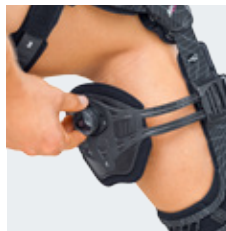

Tensioning Dial
the integrated tensioning dial enables the protection pad to be adjusted exactly to suit the individual patient.

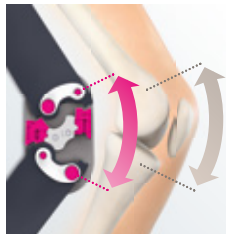


The physioglide hinge
mimics the rolling / gliding movement of the knee and ensures a secure fit.



Pre-shaped tibia pad
the winged shape of the tibia cushion moulds itself perfectly to the anatomical shape of the shin bone.

 black

More information to the exercises, lectures with international knee experts recorded at the ESSKA 2016 and much more you find on our website: www.medi.de/en

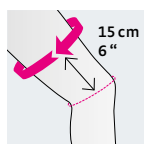

circumference thigh (cm)	31 – 37	37 – 45	45 – 51	51 – 57	57 – 65	65 – 74
size	XS	S	M	L	XL	XXL
article number	G.081◇01	G.081◇02	G.081◇03	G.081◇04	G.081◇05	G.081◇06

◇ = fill in for orientation: 2 = left, 3 = right

Following wedges are pre-installed on the brace at the point of delivery: extension 0° / flexion 90°