

medi SAS[®] 15

Shoulder abduction splint –
15° abduction



Now with breathable lower arm bag


Indications

- reconstruction after rotator cuff ruptures
- impingement syndrome
- after removal of calcium deposit
- conservative or post-operative treatment after subacromial bursitis
- conservative immobilisation of shoulder luxations / dislocations
- shoulder prosthesis insertion

Mode of action


- lateral positioning with anatomically shaped cushion to relieve the ventral limbus and the capsule
- 15° abduction position to avoid adhesions to axillary recess (frozen shoulder)
- immediate post-operative, active and passive mobilisation of hand, wrist and elbow joint by opening of forearm sling without taking off the whole support
- early functional follow-up treatment; easier post-operative physiotherapy

Product benefits

 grey




Breathable lower arm bag for improved comfort in wearing.

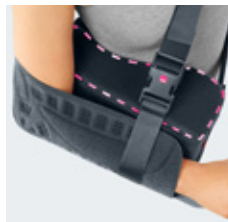

Safe and comfortable positioning by anatomically preformed pillow.



Hand rest supports the hand and allows light forearm exercises to be carried out.

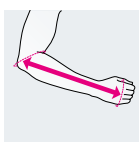

| | | | |
|-------------------------------|-----------|-----------|-----------|
| length of forearm (cm) | up to 34 | 35 – 39 | over 39 |
| size | S | M | L |
| article number | R.011.012 | R.011.013 | R.011.014 |